

High Speed Microcentrifuge















Combi-Rotor™, accepts tubes and PCR strips







ts tubes on the second of the Remaining 10:00












High Speed Microcentrifuge



- Full Color Touch Screen
- mySpin™ 1 button program recall
- COMBI-Rotor™ accepts tubes and strips
- Compact design, 9 x 12 x 7.75 in.
- Unique air-flow system, keeps samples cool



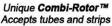



Full Color TouchScreen







Benchmark's MC-24 Touch is the first microcentrifuge in its class to feature touch screen control. In addition to the next generation control panel, the centrifuge also offers high speed (up to 13,500rpm/16,800xg) and a unique COMBI-Rotor™ that accepts both standard microtubes (1.5/2.0ml) and PCR tubes/strips (0.2ml), making this the perfect partner for any molecular biology laboratory.

The large, color touch screen shows actual speed (in rpm or rcf) along with set and remaining time. Up to 9 speed/time profiles can be stored in memory and the home screen displays a "mySpin"icon for instant recalls a designated favorite profile.

Supplied with Benchmark's unique COMBI-Rotor™, the MC-24 Touch will accept 24x1.5/2.0ml tubes, 2 strips of 0.2ml tubes (16x0.2ml), and up to 12x spin column kits without requiring separate adapters. There is also an optional 8 place rotor for 5ml tubes and cryovials. Samples are kept cool, even during extended runs at the maximum speed setting by a proprietary air flow system that moves air evenly and quietly throughout the chamber.

The MC-24 Touch is cold room safe and carries a 2 year warranty.



Technical Data:

200 to 13,500rpm 100 to 16,800 xg

Capacity: 24 x 1.5/2.0ml 2 x PCR Strips (16 x 0.2ml)

12 x Spin column

100 rpm / 100 rcf 30 sec. to 99 min. / cont. Timer:

23 x 30.5 x 19.7 cm

14 lbs / 6.4 kg Weight:

115, 60Hz or 230V, 50Hz

Ordering Information:

C2417* MC-24 Touch Microcentrifuge with 24 place COMBI-Rotor

C2417-R5 Optional 8 x 5ml rotor

C2417-R5-ADP2 Cryovial adapters for 5ml rotor, pack of 8

C1008-A5-6 Optional adapters for 0.5ml tubes, pack of 6

*115V with US plug. To order in 230V, add (-E) to the item number.

Benchmark Scientific Inc. PO Box 709 Edison, NJ 08818 Phone: 908-769-5555 Fax: 732-313-7007 Web: www.BenchmarkScientific.com Email: Info@BenchmarkScientific.com